



#### **About Us**

Aktif Dis Ticaret, was established in 1993 by the purpose of serving the Medical Sector.

Our first commercial activity was importing the document microfilm recorder -which lets archiving easy- to Turkish market.

Started to representate world-famous radiology accessories manufacturers in Turkey in 1995. Meanwhile we increased our speciality about protective lead aprons by importing and distrubiting all over to Turkey. By the experience on this area Aktif Dis Ticaret started to manufacture protective lead aprons by his own brand AKTİF® X-RAY. AKTİF® X-RAY Lead Aprons and other protective materials based on lead and leadfree are manufactured by %100 Turkish capital in Turkey and exported to world by warranty of Aktif Dis Ticaret.

We remain at your disposal by the policy of %100 Customer Satisfaction and Total Quality. AKTİF® X-RAY is registered to Aktif Dis Ticaret by Turkish Patent Institute.

Whenever you wish to get more information about our products, you can contact with us. Our team of specialists will give you the best assistance for your needs.

> Best regards, Hayri Gokturk General Manager

Tel: +90 216 474 46 48 | Fax: +90 216 474 48 19

www.kursunonluk.com.tr | info@aktifdisticaret.com





### **Coat Apron with Velcro**

This more preferred model is designed for ones who are exposed to x-ray for long term. This is an ideal model for protection against to scattered X-ray expose. It gives maximum comfort, flexibility and lightness to operators by easy wear and take-off.

Thanks to its speacial design, perfectly covering the shoulder and body and adjustable to person for maximum protecting even on shoulder bones.

#### Lead

LightWeight Lead LeadFree (optional)

### **Lead Equivalent**

Front: 0.35mm Pb or 0.50mm Pb

Size

**SMALL** 

Shoulder

Width

40 cm

**Hip Size** 

60 cm

**MEDIUM** 

Shoulder

Width

48 cm

**Hip Size** 

60 cm

**LARGE** 

Shoulder

Width

53 cm

**Hip Size** 

75 cm

(Pb)

0,35

mm

mm

mm

0,50

mm

0,35

mm

0,50

mm

Back: 0,25mm Pb



# **Coat Apron For Anesthesia**

It gives very comportable usage by its 16cm elastic belt and large velcro holders on back protection extention. You can keep it weared during your working with full of easy.

#### Lead

LightWeight Lead LeadFree (optional)

**Lead Equivalent** Front: 0,35mm Pb or 0,50mm Pb

Back: 0,25mm Pb



Size	Front (Pb)	Back (Pb)	Length	Article	
CMALL			90 cm	153	
SMALL	0,35	0,25	100 cm	154	
Shoulder Width	-	mm mm		155	
40 cm			120 cm	156	
40 Cm			130 cm	157	
			90 cm	158	
Hip Size	0,50	0,25	100 cm	159	
60 cm	mm	mm	110 cm	160	
oo ciii			120 cm	161	
			130 cm	162	
MEDIUM			90 cm	163	
Shoulder	0,35	0,25	100 cm	164	
Width	mm	mm	110 cm	165	
48 cm		"""	120 cm	166	
40 (111			130 cm	167	
			90 cm	168	
Hip Size	0,50	0,25	100 cm	169	
60 cm	mm	mm	110 cm	170	
00 (111			120 cm	171	
			130 cm	172	
LARGE			100 cm	174	
Shoulder	0,35	0,25	110 cm	175	
Width	mm		mm	120 cm	176
53 cm			130 cm	177	
) JJ (111			100 cm	179	
Hip Size	0,50	0,25	110 cm	180	
75 cm	mm	mm	120 cm	181	
/5 (111			130 cm	182	



# **Coat Apron with Belt**

This model is mounted to extention of back protection and it is ideal to distribute the load on shoulders for long-term use.

#### Lead

LightWeight Lead LeadFree (optional)

**Lead Equivalent** Front: 0,35mm Pb or 0,50mm Pb Back: 0,25mm Pb



Size	Front (Pb)	Back (Pb)	Length	Article
CMALL			90 cm	102
SMALL Shoulder	0,35	0,25	100 cm	103
Width	mm		110 cm	104
40 cm		mm mm	120 cm	600
40 Cm			130 cm	601
			90 cm	106
Him Cina	0,50	0,25	100 cm	107
Hip Size 60 cm	mm	mm	110 cm	108
ou cili			120 cm	602
			130 cm	603
MEDIUM				109
Shoulder	0.35	0,35 0,25 mm mm	100 cm	110
Width			110 cm	111
48 cm			120 cm	112
46 CIII			130 cm	113
			90 cm	114
Hip Size	0,50	0,25 mm	100 cm	115
60 cm	mm		110 cm	116
oo ciii			120 cm	117
			130 cm	118
LARGE			100 cm	119
Shoulder	0,35	0,25	110 cm	120
Width		mm	120 cm	121
53 cm			130 cm	122
33 CIII	0.50	0.35	100 cm	123
Hip Size	0,50	0,25	110 cm	124
75 cm	mm	mm	120 cm	125
/3 CIII			130 cm	126



# **Front Apron**

This more preferred model is designed for ones who are exposed to x-ray for long term. This is an ideal model for x-rays coming from front side of operators. It gives maximum comfort, flexibility and lightness to operators by easy wear and take-off.

It can be produced with belt, thick banded and with velcro upon customer request.

#### Lead

LightWeight Lead LeadFree (optional)

**Lead Equivalent** 

Front: 0,35mm Pb or 0,50mm Pb

Size	Front (Pb)	Length	Article
CMALL		90 cm	350
SMALL	0,35	100 cm	351
Shoulder Width	mm	110 cm	352
40 cm		120 cm	353
40 Cm		130 cm	354
		90 cm	355
Uin Cina	0,50	100 cm	356
Hip Size 60 cm	mm	110 cm	357
ou cm		120 cm	358
		130 cm	359
MEDIUM		90 cm	360
Shoulder	0,35	100 cm	361
Width	mm	110 cm	362
48 cm		120 cm	363
46 CIII		130 cm	364
		90 cm	365
Hip Size	0,50	100 cm	366
60 cm	mm	110 cm	367
00 (111		120 cm	368
		130 cm	369
LARGE		100 cm	371
Shoulder	0,35	110 cm	372
Width	mm	120 cm	373
53 cm		130 cm	374
33 CIII		100 cm	376
Hip Size	0,50	110 cm	377
75 cm	mm	120 cm	378
/3(111		130 cm	379



# **Double-Sided Aprons**

It gives maximum protection and preferred by operator of Kardiology section and fields where there is too much x-ray. It has 2 model with belt and witout belt giving maximum protection.

It can be weared like a jacket and its cruvaze style overlaps to eachother in order to give full protection on front side of human body and gonad region. The shoulder wadding decreases the load over the shoulders.

Additional snap festener or rubber belt is made upon request. In your orders for double-sided apron models, please indicate your preferences about snap fastener or rubber belt.

#### **Lead** LightWeight Lead LeadFree (optional)

**Lead Equivalent** 

Front: 0,35mm Pb or 0,50mm Pb

Back: 0,25mm Pb

Size	Front (Pb)	Back (Pb)	Length	Velcro	Belt	Snap
SMALL			90 cm	228	253	283
Shoulder	0,35	0,25 mm	100 cm	229	254	200
Width	mm		110 cm	230	255	201
40 cm			120 cm	612	256	202
			130 cm	613	257	203
			90 cm	232	258	204
Hip Size	0,50	0,25	100 cm	233	259	205
90-110	mm	mm	110 cm	234	260	206
cm			120 cm	614	261	207
			130 cm	615	262	208
			90 cm	235	263	209
MEDIUM Shoulder	0,35	0,25	100 cm	236	264	210
Width			110 cm	237	265	211
48 cm			120 cm	238	266	212
			130 cm	239	267	213
			90 cm	240	268	214
Hip Size	0,50	0,25	100 cm	241	269	215
110-125	mm	mm	110 cm	242	270	216
cm			120 cm	243	271	217
			130 cm	244	272	218
			100 cm	245	274	220
LARGE	0,35	0,25	110 cm	246	275	221
Shoulder Width	mm	mm	120 cm	247	276	222
53 cm			130 cm	248	277	223
			100 cm	249	279	225
Hip Size	0,50	0,25	110 cm	250	280	226
125-140	mm	mm	120 cm	251	281	227
cm			130 cm	252	282	284









Double Sided Apron with Belt

Double Sided Apron with Snap

# **Vest-Skirt Aprons**

This so chick skirt and vest combination model distributes the weight in between shoulders and waist giving you ease movement capability.

Skirt and vest model can be produced besides of our standart body size givin above according to exact body size of operators upon request.

**Lead** LightWeight Lead LeadFree (optional)

## **Lead Equivalent**

Front: 0,25 mm Pb or 0,35 mm Pb (At the overlap place in front side the lead value is double) Back: 0,25 mm Pb







Size	Front (Pb)	Back (Pb)	Vest Length	Skirt Length	Article
CNANII	0,25 mm	0.25	60 2122	Γ.Γ. avaa	415
SMALL	0,35 mm	0,25 mm	60 cm	55 cm	401
MEDIUM	0,25 mm	0,25 mm	60 cm	60 cm	416
MEDIOM	0,35 mm	0,23 111111	60 CIII	00 CIII	402
LARGE	0,25 mm	0,25 mm	60 cm	65 cm	417
LANGE	0,35 mm	0,23 111111	00 (111	05 (111	403

# **Vest-Skirt Aprons Wrap Model**

This so chick skirt and vest combination model distributes the weigth in between shoulders and waist giving you ease movement capability.

Skirt and vest overlaps on front side giving complete protection. If it is produced by 0,25mm Pb material, it becomes very light weight and gives protection of 0,50mm Pb on overlap area in front.


Skirt and West model can be produced besides of our standart body size givin above according to exact body size of operators upon request.

#### Lead

LightWeight Lead LeadFree (optional)

**Lead Equivalent** 

Front: 0,35mm Pb or 0,50mm Pb (At the overlap place in front side the lead value is double)

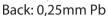



Size	Front (Pb)	Back (Pb)	Vest Length	Skirt Length	Article
CNAALL	0,35 mm	0.25	60		430
SMALL	0,50 mm	0,25 mm	60 cm	55 cm	420
MEDIUM	0,35 mm	0.25 mm	60 cm	60 cm	431
MEDIOM	0,50 mm	0,25 mm	00 CIII	o citi   oo citi	421
LARGE	0,35 mm	0,25 mm	60 cm	65 cm	432
LANGE	0,50 mm	0,2311111	00 CIII	05 (111	422



### **Double-Sided Apron Wrap Model**

This wrap model apron gives perfect protection on front side of the body because it covers front side of the body completely. For the two sided x-ray expose, overall protection level is keeped as 0,25mm Pb at single front part, so it gives 0,50 mm Pb protection for long term use by a lighter weight apron.

It gives maximum protection and preferred by operator of Cardiology department and fields where there is too much x-ray. It has 2 model with belt and witout belt giving maximum protection.

It can be weared like a jacket and its cruvaze style overlaps to eachother in order to give full protection on front side of human body and gonad region. The shoulder wadding decreases the load over the shoulders.

#### **Lead** LightWeight Lead LeadFree (optional)

**Lead Equivalent** 

Front: 0,35mm Pb or 0,50mm Pb

Back: 0,25mm Pb



Size	Front (Pb)	Back (Pb)	Length	Article
CMALL			90 cm	520
SMALL Shoulder	0.35	0,35 0,25 mm mm		521
Width	mm			522
40 cm		'	120 cm	523
40 CIII			130 cm	524
Hip Size			90 cm	525
90-110	0,50	0,25	100 cm	526
cm	mm	mm	110 cm	527
			120 cm	528
			130 cm	529
			90 cm	530
MEDIUM	0.35	0,35 0,25 mm mm	100 cm	531
Shoulder	1		110 cm	532
Width			120 cm	533
48 cm			130 cm	534
Hip Size			90 cm	535
110-125	0,50	0,25	100 cm	536
cm	mm	mm	110 cm	537
		'''''	120 cm	538
			130 cm	539
LARGE			100 cm	541
Shoulder	0,35	0,25	110 cm	542
Width	mm		120 cm	543
48 cm			130 cm	544
			100 cm	546
Hip Size	0,50	0,25	110 cm	547
125-140	mm	mm	120 cm	548
cm			130 cm	549



# **Dental Aprons**

Dentists uses medical radiological examination technics as Panoramic, Periapical, Sefolometric, Sinus graphics very often for diagnonistic. That's why it is necessary to use lead apron and thyroid collars during dental treatment. Speacialy for children and pregnants it is not allowed to make shooting without lead apron.

#### **Lead** LightWeight Lead LeadFree (optional)

**Lead Equivalent** 

Front: 0,35 mm Pb or 0,50 mm Pb

Back: 0.25 mm Pb

Size	Front (Pb)	Back (Pb)	Length Width	Article
CAAAII	0,35 mm	0.25	60,460	341
SMALL	0,50 mm	0,25 mm	60 x 60	320
MEDIUM	0,35 mm	0.25 mm	70 x 60	342
MEDIOM	0,50 mm	0,25 mm	70 X 60	321
LARGE	0,35 mm	- 0,25 mm 75 x 75	343	
LANGE	0,50 mm		/5 X /5	322





# **Patient Protection Apron**

These model is intentionally produced for patient under dental treatment. This model is ease of use and has thyroid protector included giving maximum protection to user.

**Lead**LightWeight Lead
LeadFree (optional)

**Lead Equivalent** Front: 0,35 mm Pb or 0,50 mm Pb

Size		Dimensions	Article
STANDARD	0,35 mm	60 x 100 cm	323
STANDARD	0,50 mm	00 X 100 CIII	340

# **Surgical Aprons**

It is weared under the surgical gown giving protection to operator in case of continius shooting. Doctors can not take off by themselves the protection apron because they can not take a risk of hand sterilisation. That's why an assistant can take off the velcro band, then the protection lead apron drop down. No need to use this apron if the shooting is finished.

**Lead** LightWeight Lead LeadFree (optional)

**Lead Equivalent** 0,25mm Pb or 0.35mm Pb





Size	Front (Pb)	Article
50 x 80 cm		405
60 x 90 cm	0,25	406
60 x 100 cm		407
60 x 110 cm		408
60 x 120 cm		409
50 x 80 cm		410
60 x 90 cm	0,35	411
60 x 100 cm		412
60 x 110 cm		413
60 x 120 cm		414

# **Thyroid Collars**

It is produced in order to protect the thyroid region which is one of the most effected part in human body from x-ray. It has 5 models: Standard, Hat Model, Large Hat Model, Panoramic and Children.

It is recommended operators to have their own private thyroid collar as a matter of hygiene (sweating).

The standart model can be mounted to double-sided aprons by snap fastener. It should be sewed to only front side closed aprons.

#### **Lead** LightWeight Lead LeadFree (optional)

**Lead Equivalent** 0,50 mm Pb

Model	Pb/mm	Article
Panoramic		316
Large Hat Style		317
Standard	0,50 mm	318
Hat Style		319
For Children	]	328



**Large Hat Style** 



**Standard** 



**Hat Style** 

# **Aprons for Children**

Children aprons are produced mostly as front and back protected. They are produced with belt fitting to body well and gives maximum safety and easy of use.

#### **Coat Apron for Children**

### **Lead** LightWeight Lead LeadFree (optional)

**Lead Equivalent** Front: 0,50 mm Pb Back: 0,25 mm Pb

Size	Front (Pb)	Back (Pb)	Article
3-5 YEARS		0.25	501
6-8 YEARS	0,50 mm		502
9-12 YEARS		0,25 mm   	503
SET			504





# **Double-Sided Apron for Children**

These models are ease to wear and take off. They are tightened by velcro simply. They gives maximum protection and safety for children.

LeadLead EquivalentLightWeight LeadFront: 0,50 mm PbLeadFree (optional)Back: 0,25 mm Pb

Size	Front (Pb)	Back (Pb)	Article
3-5 YEARS	0,50 mm		505
6-8 YEARS		0.25	506
9-12 YEARS		U,25 MM	507
SET			508

# **Multi purpose Protective Covers**

These are ease of use and can be tightened to weist by means of velcro belt. They can be used for gonad protector or many other purpose.

**Lead** LightWeight Lead LeadFree (optional)

**Lead Equivalent** 0,5 mm Pb

Size	Pb	Article	
20x25 cm	0,50 mm	301	
30x30 cm		302	
37x40 cm		303	
40x50 cm		304	
SET		305	





#### **Protective Covers**

Any of the protective covers below can be used for all radiological application. The important point is select the proper size for your need.

# **Lead** LightWeight Lead

LeadFree (optional)

# **Lead Equivalent** 0,35mm Pb or 0,5 mm Pb

Size	Pb	Article	
60x60 cm		306	
90x60 cm	0,35 mm	307	
120x60 cm	]	308	
60x60 cm		309	
90x60 cm	0,50 mm	310	
120x60 cm	]	311	

# **Ovary Shields and Gonad Shields**

#### **Ovary Shields**

Ovary is the most sensitive regions of our bodies -where our genitals are located- to X-RAY. AKTIF X-RAY oavry shileds are ease of comfort and clean and starts from 1.00 mm Pb lead level providing the maximum protection. It is produced as 3 different models for 0-9 years, 9-16 years and over 16 years users.

#### Lead

LightWeight Lead LeadFree (optional)

# **Lead Equivalent** 1,0 mm Pb

Size	Pb	Article
Small	1,00 mm	312
Medium		313
Large		314
Set		315



#### **Gonad Shields for Males**

It must be used to protect gonad region -where our genitals are located- to X-RAY.

It is produced as 3 different models for 0-9 years, 9-16 years and over 16 years users.

#### Lead

LightWeight Lead LeadFree (optional)

# **Lead Equivalent** 1,0 mm Pb

Size	Pb	Article
Small	1,00 mm	324
Medium		325
Large		326
Set		327

# **Leg Protectors**

Leg protector is produced according to customer request and it has velcro fastener closer that enable to use in every size.



#### Lead

LightWeight Lead LeadFree (optional)

#### **Lead Equivalent**

0,50 mm Pb



# **Special Models**

We can produce high quality and perfectly produced special products for patient and operators.

#### **Arm Protectors**

We can produce arm protectors upon request. 6 different colors and lead levels can be selected.



**Arm Protectors** 

#### **Back Relief Belt**

Lead aprons are heavy protectors and belts are loosing its function in time, that's why Back Relief Belts are needed. These special design belts can save its flexibility for a long time but its similar products loses the flexibility in a very short time. Wearing and taking off is very easy. It is adjustable and safety product with velcro and claps. They can be produced as 4 different sizes (S,M,L,XL) and in different colors upon customer request.

Features: S,M,L,XL. Claps and Velcro.



#### **Breast Shield**

It is used for spinal cord shots. Height can be adjustable with velcro band. It does not contain lead on neck part for full spinal cord shots. It is ideal for scollosis works.



**Lead Caps** 

#### **Racks**

It is important to save aprons in proper condition in order to last their life. At the same time, easy reach, user friendlyness, small size and easy mounting are essential points of saving style.

#### Wall Mounted Apron Rack: (2 arms - 3 arms - 5 arms)

This user friendly and ergonomic hanger help to long last your lead apron's life. Special hooks can be mounted to this model for skirts.



#### **Mobile Apron Rack**

Capacity is 5-10 aprons with Natural Wooden Hangers. Please indicate the quantity of Natural Wooden Hangers with your orders.



Model	Article
Wall Mounted Rack For Gloves	400008
Wall Mounted Rack For Two Aprons	400025
Wall Mounted Rack For Three Aprons	400105
Wall Mounted Rack For Five Aprons	400104
Mobile Rack For Five Aprons Stainless Steel	400011
Mobile Rack For Ten Aprons Stainless Steel	400013
Natural Wooden Hanger	400015
Mobile Apron Rack	400111
Lux Mobile Rack with Shelf	400112



# **Natural Wooden Hanger**

Strong wooden body, stainless steel hook part and proper design long lasts your apron's life. Capasity is only one apron. Color is teak.



#### **Lux Mobile Rack with Shelf**

This mobile model has metal foots and strong wheels making it to move easly. 10-12 aprons can be hanged. It has multi purpose use shelves. Please indicate the quantity of Natural Wooden Hanger with your orders.



# **Lead Surgical Gloves**

This gloves are made of natural dry rubber mixed with lead oxide. The mix of these two raw materials grants high elasticity and softness combined with the best x-ray protective curve without negative peaks tipical of other protective materials. Each single gloves, bisedes, is controlled and tested by specialized operators that verify the full conformity to the Standards requisities and to the strictly requirements defined by our Quality Management System.

#### **Models:**

XP-1\* Light Type XP-2\*\* Medium Type XP-3\*\*\* Heavy Type

#### Where to use:

**XP-1\* Light Type:** Gives maximum touch sense, it is self-lubricated (ease wearing), used more at angiography, urographi, hemotology, gastrology, light-simple surgical operations.

**XP-2\*\* Medium Type:** Standart thickness, self-lubricated (ease wearing), used more at urology, oncology centers, dentistry, surgery and cardiology.

XP-3\*\*\* Heavy Type: It has thicker materials and proper for traumatology and nuclear field.

Size	XP-1* Light	XP-2** Medium	XP-3*** Heavy
6 ½	750	756	762
7	751	757	763
7 ½	752	758	764
8	753	759	765
8 1/2	754	760	766
9	755	761	767

Absorbtion and Screening Percentage	60 Kv	80 Kv	100 Kv
XP-1* Light	36%	31%	27%
XP-2** Medium	61%	54%	47%
XP-3*** Heavy	72%	65%	59%





2 stars

3 stars

#### **Lead Gloves**

Aktif X-ray Lead Gloves are flexible, homogenius and gives maximum safety to operator with 0,50mm Pb lead value.

Size	Pb	Article
Standard	0,50 mm	772



# **Lead Eyewear**

Our lead eyewear is well-selected and gives the highest protection and comfort of use to operators. Panaromic type eyewear is adviced to be used in areas where there is radiation expose risk from sides of the head. PS: Prescription panoramic protection eyewears can be produced.



Model	Pb	Article
Panorama	0,75 mm	700
Ultralite 99	0,75 mm	701
Fitover	0,75 mm	702
Standard	1,00 mm	703
Ultralite 9935	0,75 mm	704
Ultralite 9941	0,75 mm	705
Alumilities	0,75 mm	706
99AL	0,7,5 111111	, 30

#### **Panoramic Half Protection Mask**

A Mask of leaded acrylic protection -covering over to almost complete surface of the face- can be applied.

It is ideal for minumum bad-effect of radiation.

It can be cleaned very easly.

It must be keeped away from high pressure and temperature.

#### **Panoramic Complete Protection Mask**

Leaded acrylic mask protects the coplete face from forehead to chin.

It has 20cm length.

Adjustable headband gives use of comfort. It must be keeped away from high pressure and temperature.





**Panoramic Half Protection Mask** 



**Panoramic Complete Protection Mask** 

### **Multi Purpose Lead Curtains**

Lead curtain is produced in order to protect humans against scattered x-rays. They can be used and applied with ease.

**Applications:** Aisle, cabin entry, seperating the working areas, covering windows and many other purpose and places. Mobile lead curtain can be slided to different ways with ease giving you usability of multi-purpose positions.

#### Lead

LightWeight Lead LeadFree (optional)

#### **Lead Equivalent**

0,25 mm Pb, 0,35 mm Pb, 0,50 mm Pb

#### **Curtains for Table**

Standard 12x80 size protective curtains can be opened and closed by overlapping each part to others system. In operation rooms it can be applied to C armed rontgen tube. It can be produced in requested sizes. It can be produced with 0,25/0,35/0,50 lead equivalent. Select your color.

Size	Pb
12x80 cm	0,25 mm 0,35 mm 0,50 mm
Other size on demand	





#### **Lower Body Shields for Operation Tables**

It can be mounted to operation table's accessories rail and protecting operator agains secondary radiation coming through under secondary table.

Table mounted model is designed to protect doctor and other stuff against to radiaton during operation. Protective lead material's size is 116 cm x 76 cm, lead equivalent is 0,50 mm Pb. Metal part is not magnetic. It can be mounted to both sides of operation table. It can move to front-back direction with ease on table accessory rail giving maximum comfort and safety.

#### **Technical Features**

**Lead Equivalent:** 0,5 mm Pb flexible curtain **Dimensions: Length:**116 cm x Heigth:76 cm

**Color:** Grey

Metal Part: Eloxal Covered Aluminum



#### **Standard Paravane**

Dimension of standart lead paravane is 100cm x 190cm. Strong wheel system and eye height level lead glass give relax movement capability and wide viev perspective to operator.

Model	Size	Pb	Glass Size	Article
Standard Paravan	95 x 195 cm	2,0 mm	20 x 20 cm	720



# **Standard Paravane with Large Glass**

Although the technical characteristics and sizes same as the standard paravane with a larger glass offers the user a wider viewing angle.

Model	Size	Pb	Glass Size	Article
Standard Paravane with Large Glass	95 x 195 cm	2,0 mm	50 x 80 cm	721



#### **Half Glazed Acrylic Paravane**

Diagnostic using X-RAY requires personal to be at examination room. In this case personal is exposed to radiation. Places like operation rooms, intensive care units, Cardiac Units-where X-RAY applied- needs to be protected. In this places, this half glazed paravane-which has large angle of view- usage is adviced beside protective overall. Steel parts are stainless. Acrylic glass can be asjusted upon requested height.

Model	Size	Pb	Glass Size	Article
Half Glazed Acrylic Paravane	180 x 70 cm	2,0 mm	60 x 70 cm	722



**Half Glazed Acrylic Paravane** 

# X-RAY DOES NOT JOKE!

X-Ray Protection
Coat Aprons

Coat Aprons for Children Surgeon Aprons

Vest & Skirt

**Double-Sided Aprons** 

Multipurpose Protective Covers

**Dental Aprons** 

**Thyroid Collars** 

**Protective Cover for Patients** 

Ovary Shields

Gonad Shields

Lead Eyewears

Lead Gloves

Lead Gloves for Surgery

Racks for Lead Aprons

Lead Markers

**Intensifying Screens** 

Compensating Filters

**ID Cameras** 

Digital ID Cameras

Cabinets and Tables for X-Ray Department

**Quality Control Equipments** 

Negatoscopes

X-Ray Films

X-Ray Cassettes

Mammography Cassettes

Multiformat CRT Cassettes

Carbonfibre Cassettes

Radiotheraphy Cassettes

Cassette Divider 1.0 Pb

Flexible Cassettes

**Curved Cassettes** 

Cassette Holders

Grids

Pneumocolon for Double Contrast

**Positioning Aids** 

Patient Table Pads

Elusson Ultrasound Gel

Screen Cleaners
Hilus Filters

Kusbakisi Cad. Asuroglu Sitesi D-Blok 29/1-2 Altunizade, 34662 Uskudar, İstanbul-Turkey Tel: +90 216 474 46 48 Fax: +90 216 474 48 19 www.kursunonluk.com.tr